

Flow Cytometry: Principles and Instrumentation

Rafael Nunez*

Associate Laboratory Member, Memorial Sloan-Kettering Cancer Center, Box 98, 1275 York Avenue, New York, NY 10021, USA

Introduction

The standard benchtop flow cytometer is very similar to a hematology cell counter. In fact, flow cytometers can trace their origins back to the early hematology counters. Unlike their earlier counterparts that used electronic impedance to measure particles in a fluid stream, today's modern flow cytometers use an illuminating light source usually consisting of a laser or arc lamp. The majority of instrument manufacturers employ an air-cooled argon gas laser emitting a monochromatic beam of light fixed at 488 nm at 15 mW of power. As particles or cells flow in single file past the intersection of the light beam, light is scattered in various directions. If there is a fluorochrome labeled monoclonal antibody associated with the cell, it becomes excited by the laser and a fluorescent emission results. The resulting signals are processed to gather information about the relative size of the cell (forward light scatter-FSC), its shape or internal complexity (side light scatter-SSC) as well as a diversity of cellular structures and antigens (fluorescence).

The cytometer itself is set up and monitored routinely with a quality control program utilizing a series of unlabeled and fluorescently labeled calibration particles as a reference check on instrument alignment and sensitivity performance. Therefore, every cytometric procedure performed is a reflection of the quality of the instrument, sample preparation, and the operator, and ultimately reflects on the institution.

FACSCalibur™

The FACSCalibur uses an air-cooled argon gas laser regulated at 15 mW of power and a fixed wavelength emission of 488 nm. The instrument is capable of detecting six parameters: Forward Scatter (FSC), side scatter (SSC) and three fluorescent emissions (green, yellow-orange and red) utilizing the first laser. A smaller diode laser emitting red light at 635 nm is used to excite compounds that fluoresce above 650 nm and are detected with the fourth fluorescent detector (photomultiplier tube). Figures 1-3 illustrate the general system design and component layout of a flow cytometer. Figure 1 illustrates the fluidics of a cytometer. Figure 2 shows the optics of a cytometer with two lasers and the capability to detect up to six parameters. Figure 3 demonstrates the process of acquisition and

analysis of a sample and the capabilities to save the data on disks as well as to analyze the data and produce printouts and/or slides.

There is a large number of fluorogenic dyes available that can be covalently coupled to antibodies or used for intracellular staining of various cellular constituents. The most widely used fluorochrome for immunophenotyping is fluorescein isothiocyanate (FITC) and is detected in the green channel or FL1. Phycoerythrin (PE) emits in the yellow-orange region and is detected in FL2, while Propidium Iodide is a DNA specific dye that can be measured in either the orange or red channel (FL2 or FL3). Additional dyes like PerCP, TriColor and Red613 can be excited by 488 nm and emit in the red wavelengths.

FACSCComp™ software is a software application for automating the setup and electronic compensation necessary when using multiple fluorochromes. FACSCComp is used in conjunction with Calibrite™ beads which are polystyrene beads impregnated with the appropriate fluorochrome. The software informs the user when to insert the required tube, the compensation and sensitivity being performed automatically. Results are displayed on the video monitor, saved by the system, and a printout is available as a record for quality control purposes.

Unlike the traditional benchtop cell analyzer, the FACSCalibur can separate cells into different collection containers based on their physical properties. This process of selection is the basis for cell sorting (see below). When only analyzing cells, samples are consumed and discarded. The instrument is designed with a closed fluid system, thereby making it ideally suited for working with potentially biohazardous samples, *i.e.*, human blood. Working with blood samples must be authorized and the appropriate precautions implemented.

The FACSCalibur can analyze cell suspensions at the rate of several hundred cells per second. Ideally, each sample to be analyzed should contain a cellular concentration of between 10^5 - 1×10^6 . Typically, investigators acquire 5,000 to 10,000 cells per sample file; for sperm analysis up to 100,000 cells could be collected. Larger data files have greater storage requirements and generally require storage media such as Magneto-optical (MO) disks. Normally, data are first saved directly to the hard disk of the PowerPC Macintosh computer. After backing up to a 1.44 MB floppy or 230 MB rewritable optical disks, the data can be transferred to a FACS™ analysis data station. The off-line work station should contain a copy of the CellQuest™ data analysis software. Moreover, excellent public domain software for the PC (WinMDI™) and commercially available programs such as FlowJo™ for the Macintosh are available from sites on the internet.

The FACSCalibur is easy to use, and in our facility, the instrument is operated by the individual researcher who has participated and successfully completed the training course. The facility maintains a schedule sheet for reserving analysis time on a first-come basis. A block of time can be reserved in advance for sample acquisition or data analysis.

*For correspondence. Email nunezn@mskcc.org; Tel. (212) 639-6392; Fax. (917) 432-2333.

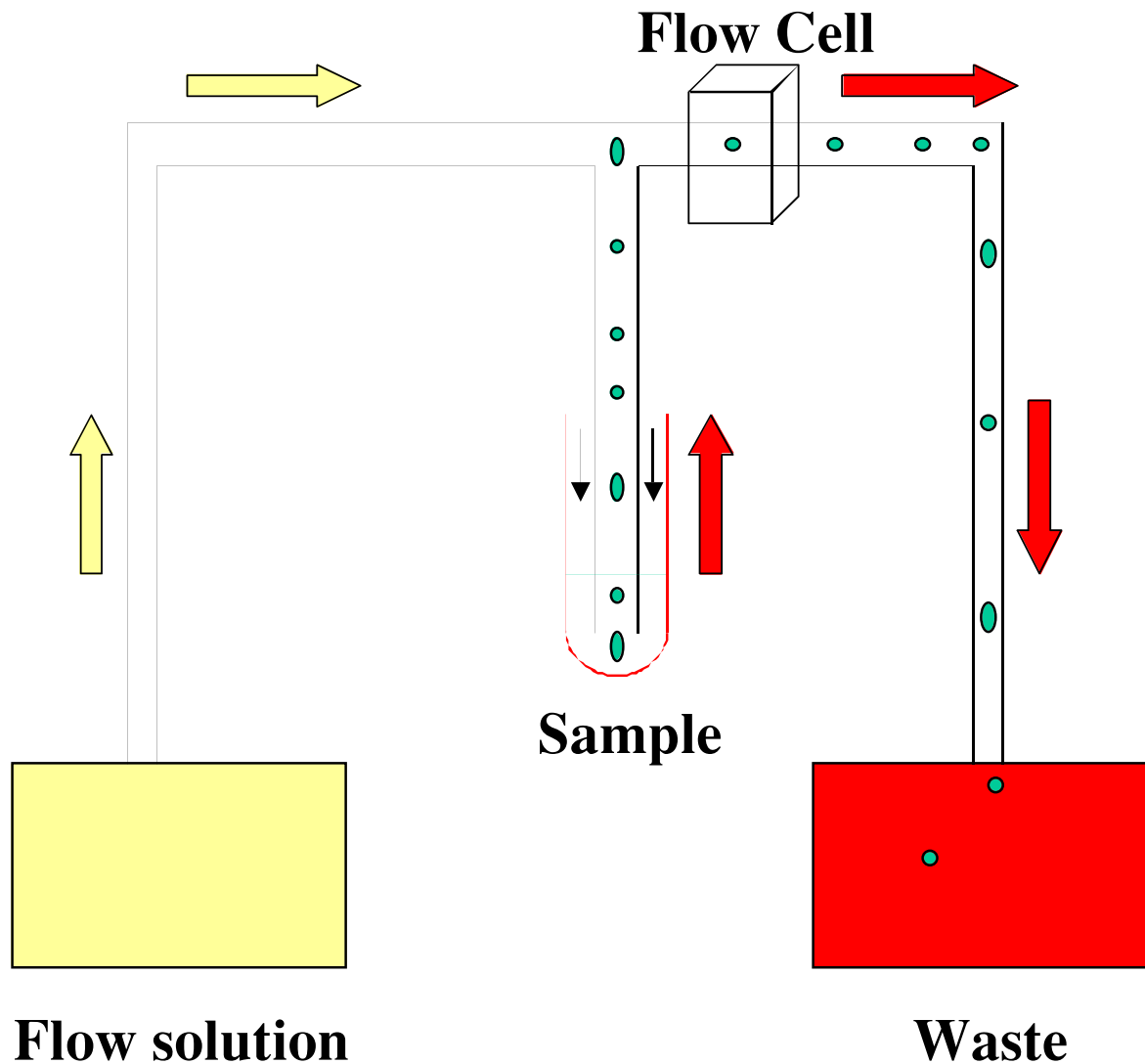


Figure 1. Fluidics. The fluidics of a cytometer.

For additional information on the FACSCalibur instrument, reagents and software, BDIS maintains a web site at: <http://www.bdfacs.com/literature/brochures/Caliburimages/CaliburBro.gif>

BD LSR

The BD LSR™ is the newest instrument in the flow cytometry line. Introduced in August 1999, this 8 parameter (6 color and two scatter measurements) is a dual laser, benchtop instrument offering simultaneous 488 nm and UV excitation. It features an air-cooled 8 mW, 325 nm He-Cd laser in addition to the 15 mW, 488 nm argon-ion laser. Building upon the ease of use of the FACSCalibur, the BD LSR simplifies routine UV applications from the more complex cell sorter instruments to a high performance, compact benchtop flow cytometer. The instrument offers

software instrument control, push-button fluidics and fine adjustment sample flow rate control. The optical system is alignment-free and permits the use of UV-excited dyes such as Hoechst 33342, DAPI and Indo-1. Our facility uses the LSR for multicolour analysis (8 parameters).

Characteristics of a Cell Sorting Instrument and Potential Applications

A cell sorting instrument is designed first and foremost as a very high-speed sorter and analyzer. It offers the sort performance needed for the quantification and/or isolation of rare cell populations with high recovery and purity. The design specification for a cell sorting instrument includes speeds at or above 25,000 cells per second, and performance approximately 50-100 times faster than existing instruments. Sorted cells must be viable with purity

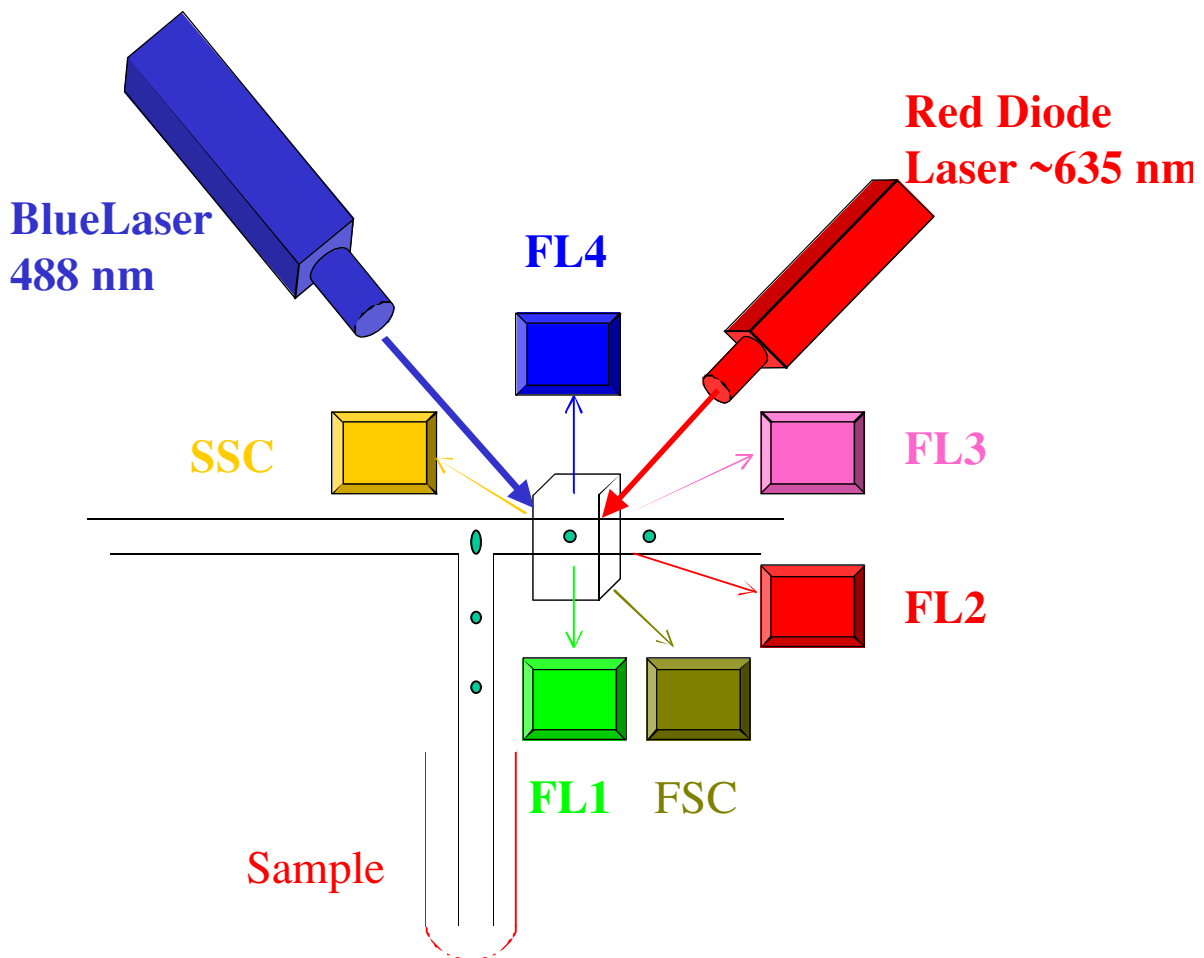


Figure 2. Optics. The optics of a cytometer equipped with two lasers and the detectors for six parameters.

and recovery better than 95 percent at any speed.

Moreover, current performance specifications for a sell sorting instrument include a speed of 50,000 cells per second, with purity and recovery greater than 99 percent at that speed.

A high speed and sensible cell sorting instrument equipped with a single-cell sorting module permits sterile cloning of single cells and allows sorting of single cells directly into tubes for subsequent PCR enzymatic amplification of DNA or RNA from sorted cells. Also on-line are special home-made, high-resolution time-of-flight sizing systems capable of sizing cells and subcellular organelles as small as 0.3 microns in length or diameter, and a multiparameter high-speed "rare-event" analysis system ("HIGHSPEED") capable of analyzing cells at rates in excess of 100,000 cells per second allowing analysis of rare cell subpopulations as small as 0.0001 percent with sample rates of more than 10^9 cells per hour.

FACSCalibur Sorting Option

As previously mentioned, the FACSCalibur is equipped with the optional sorting module. The system employs a

mechanical catcher tube assembly within the closed fluid system to divert cells to the collection tubes. This method produces no biohazardous aerosols, is ideal for pathogenic samples, but is considerably slower in sorting. For sorting of non-biohazardous samples where the sorting criteria cannot be met on the FACSCalibur, a larger multiparameter jet-in-air cell sorter should be used.

The FACSCalibur can acquire data and sort cells. While acquisition can be performed at rates over 1,000 events per second, the sorting of cells is limited to approximately 300 per second. Non hazardous living cells can be sorted, and may be recovered in gnotobiotic ("sterile") form for subsequent *in vitro* functional studies.

Although sample acquisition on the FACSCalibur is performed by the researcher, sorting of cells is a more complex procedure requiring precise alignment. Because of its open design, which permits great flexibility in configuration, its delicate components are much more easily damaged by inexperienced operators, and quite expensive to replace. Therefore operation of the FACSCalibur for sorting purposes is restricted to trained operators.

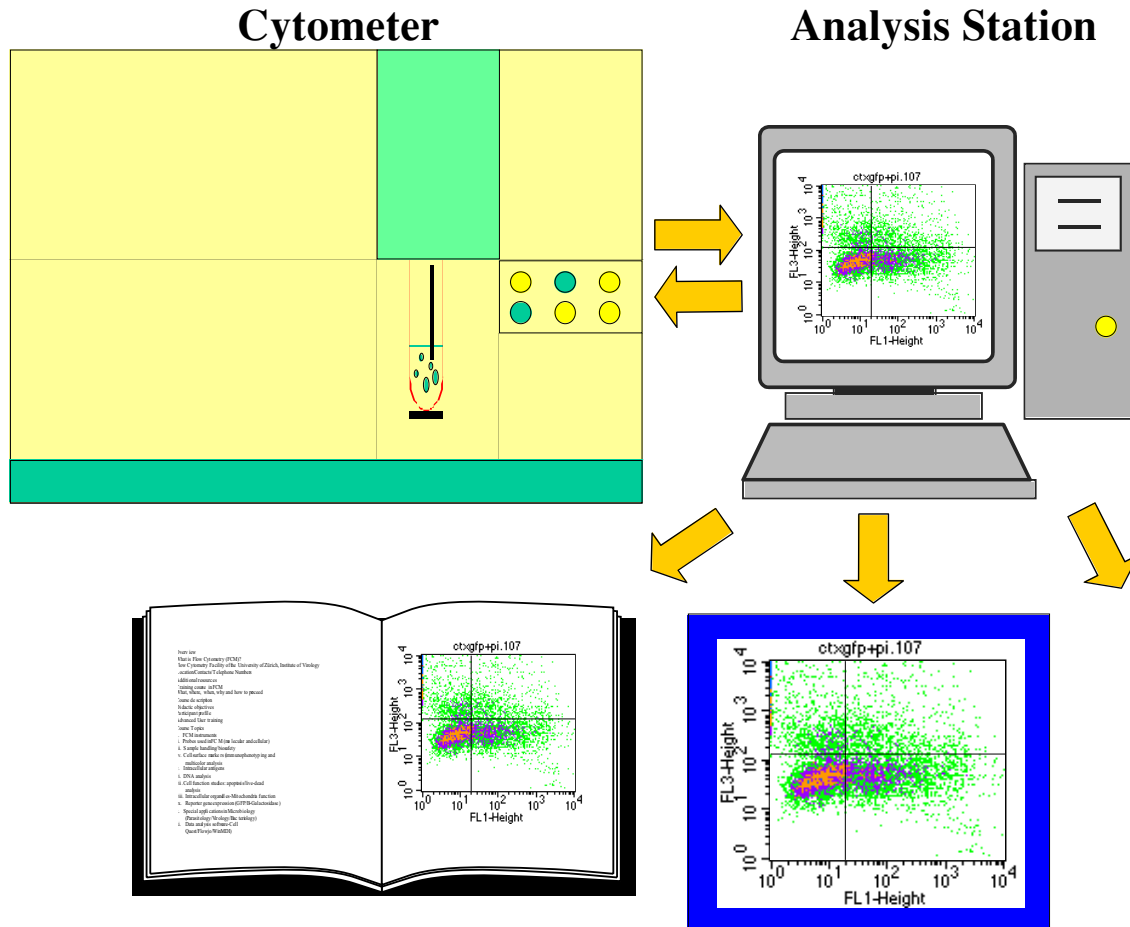


Figure 3. Acquisition and analysis. The process of acquisition and analysis of a sample and the capability to save the data on disks as well as to generate print outs and/or design slides.

FACS Vantage™ and Vantage SE

The Becton-Dickinson FACS Vantage Cell Sorter instrument is capable of sorting cells at a much higher speed than the FACSCalibur. Unlike the Calibur which has an enclosed sorting system, the FACS Vantage uses a jet-in-air system to achieve greater sort speeds. With robust analytical capabilities and flexibility in changing the optical configuration and adding of lasers, a wider range of fluorochromes can be used for analysis and sorting applications (Figure 4). Digital hardware (DIVA) was incorporated to make the sorter more user-friendly.

Beckman Coulter

Beckman Coulter manufactures the COULTER® EPICS® XL/XL-MCL FlowCytometry System for bench top operations and the “EPICS® ALTRA™ and HyPerSort™ High-Performance Cell Sorting System with capability for high speed sorting applications.

Coulter® EPICS® XL/XL-MCL FlowCytometry System

The Coulter® EPICS® XL Flow Cytometer combines the analytical power of a research cytometer into a compact, ergonomically designed clinical analyzer. The XL System features the capability to analyze up to four colors of immunofluorescence from a single air cooled laser. Other multi-color applications include: multiparametric DNA analysis, platelet studies, reticulocyte enumeration, cell biology/functional studies as well as a broad range of research applications.

For research applications, the user has the ability to change the optical filters for each fluorescent measurement channel. For high volume users, the XL-MCL System offers walk-away sample handling with the Multi Carousel Loader (MCL). This device incorporates positive bar-code identification and vortex mixing prior to sample aspiration.

The XL uses Digital Signal Processing (DSP) for reliable linearity and drift-free amplification and compensation. Multi-fluorochrome compensation is

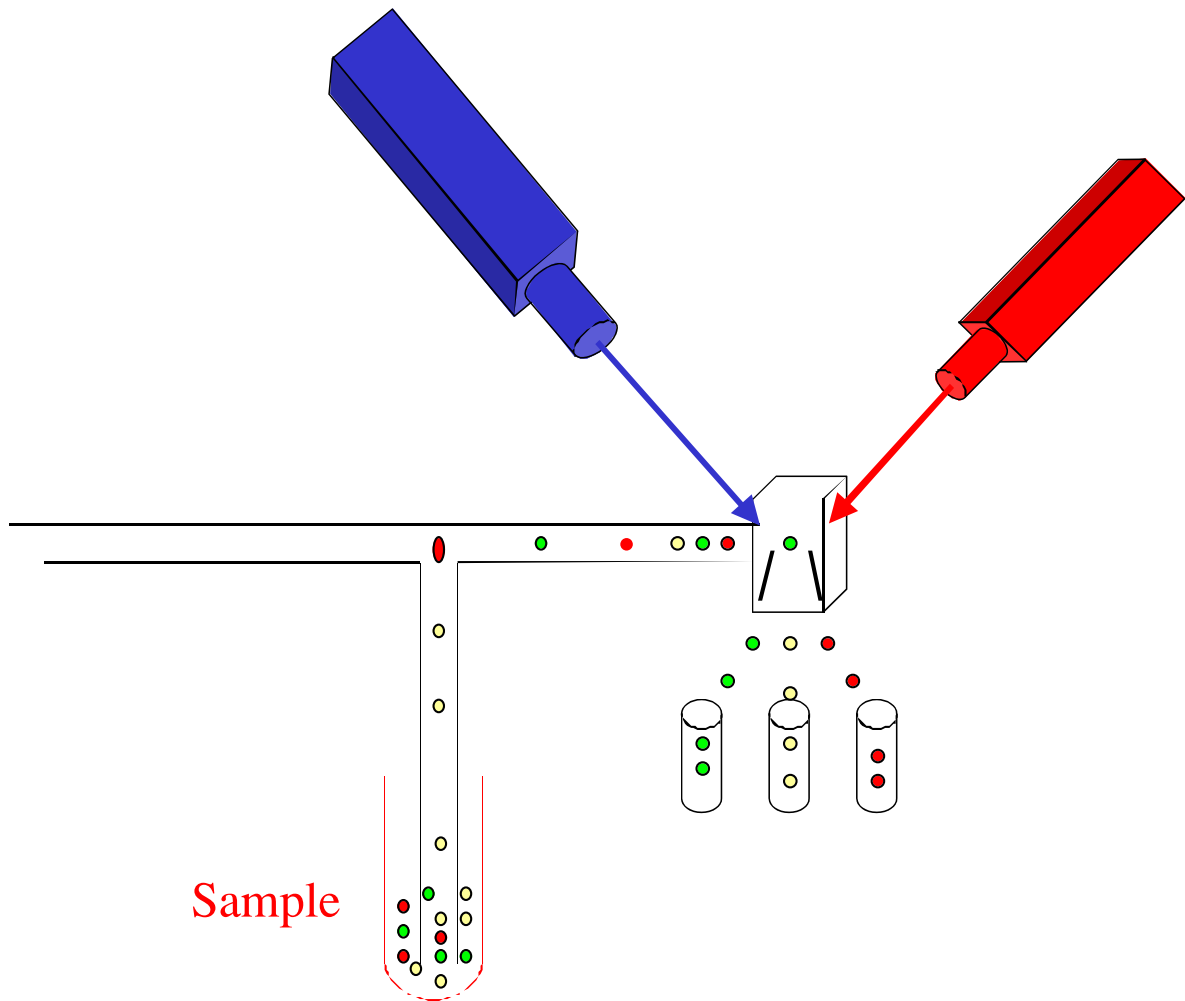


Figure 4. Sorting. This figure illustrates the general system for sorting with a flow cytometer.

achieved by a digital compensation matrix. The new XL SYSTEM II™ Software fully automates instrument set-up and compensation for 2, 3 and 4-color applications. Optional features in the software include tetraONE, the 4-color, one tube, automatic analysis algorithm and reticONE, an analysis algorithm for reticulocyte analysis.

EPICS®ALTRA™ and HyPerSort™ High-Performance Cell Sorting System

The EPICS ALTRA flow cytometer represents a new generation of cell sorters from Beckman Coulter. This high-performance instrument is a powerful tool for innovative research in immunobiology, cell physiology, molecular biology, genetics, microbiology, water quality and plant cell analysis. High fluorescence sensitivity with quartz SortSense™ flow cell design, for sense-in-quartz, jet-in-air operation. The system features the capability to analyze and sort on up to 6 colors simultaneously while performing complex multi-parameter applications. These include DNA cell-cycle analysis, quantification, functional studies, chromosome enumeration and physiological measurements.

The standard ALTRA System allows rapid separation of large numbers of specific cell populations with high purity, recovery and yield. Sorting at rates up to 10,000 cells per second (with air-cooled lasers) and up to 15,000 cells per second (with water-cooled lasers) is possible. The system features eight sort modes. Four modes optimize between purity and yield for two sorted populations, including a special enrichment mode. Three modes feature the new ALTRASort™ mixed-mode that isolates a high purity population to the left and captures the remainder of the same population to the right. This is ideal when working with rare populations and/or limited sample volume.

The AccuSort™ mode is for accurate counting when sorting individual cells for cloning with the Autoclone® sorting option or for sort matrix verification. The optional HyPerSort™ System incorporates new sorting technology that provides true high speed sorting, while maintaining sensitivity with the SortSense flow cell. HyPerSort technology uses high pressure and high frequency to maximize purity, recovery and yield at data rates up to 30,000 cells per second.

This user-configurable optical platform permits choosing from an extensive range of laser options, allowing

almost any combination of excitation wavelengths. The FlowCentre™ workstation, a high performance multi-media platform incorporates Intel® Pentium® processing chip technology. The system features EXPOTM Software for instrument control, data acquisition and analysis in a Windows 95 environment.

MoFlo® Cytometers

Cytomation <http://www.cytomation.com/index.htm> has developed diverse systems for flow cytometry. The MoFlo® MLS Flow Cytometer is a high-performance cell sorting system for high speed sorting applications. In addition, Cytomation produces a MoFlo® BTA cytometer for bench top analysis, the MoFlo® BTS benchtop sorter and the MoFlo® SX cytometer for X/Y sperm sorting separation.

The MoFlo incorporates parallel processing electronics first developed at the Lawrence Livermore National Laboratory (LLNL) and was used successfully as a chromosome sorter for the Human Genome Project. It was designed first and foremost as a high-speed sorter and analyzer, offering the sort of performance needed for the quantification and/or isolation of rare cell populations with high recovery and purity at speeds in excess of 25,000 cells per second. Current performance specifications claimed by the manufacturer are a sorting speed of 50,000 cells per second, with purity and recovery greater than 99 percent.

Cytomation offers a number of accessories for the MoFlo MLS cell sorter. 4Way™ Sorting allows researchers to sort up to 4 subsets of cells at one time, using a combination of scatter and fluorescent parameters. This new mode of sorting allows sorting of more subsets faster, and with fewer cells lost to waste.

To ensure accurate cell sorting Cytomation offers Sortmaster™ an automated system that determines the droplet delay for each sort. The system monitors the droplet break-off during the run, ready to interrupt the sort and determine a new droplet delay if an instability occurs. If a stable break-off point is found, the sort resumes. If not, the sort will be aborted and the operator notified by an audible alarm and/or pager.

CyCLONE® is an automated high-speed cloning device that sorts into trays of 96, 384 or 1536-well configuration. In addition, the software that controls the CyCLONE unit can be configured for any user specified pattern. A 96-well tray can be filled with single cells in under 50 seconds.

Specifications and Characteristics

Purity, recovery and yield of sorted cells are important parameters that researchers need to be concerned with. The ultimate goal is to obtain the highest values for each of these parameters. In the real world, there is usually a trade-off between these criteria, especially when increasing the sort speed. The MoFlo MLS is capable of generating outstanding purity results such as these at 25,000 events/second. MoFlo's electronics technology increased speed of sorting without compromising recovery and yield, making the customary problem of turning up the sort speed only to degrade recovery and yield a problem of the past.

One measure of the precision of a cell sorter's

performance is the coefficient of variation, simply abbreviated as CV. This is the standard deviation of the mean over the mean. It is a function of how well the optical elements are aligned with the sample stream. When multiple lasers are focused through a single optical element, only one beam can traverse the optical center of the lens. The other beam geometries may be compromised. This situation is best seen by noting the larger CVs generated by the off-axis beam vs. the on-axis beam. The MoFlo eliminates this situation by using separate optical elements for each beam. All beams are spatially separated and optically optimized. New software (Summit®) improved the analysis capacity of the MoFlo equipment.

The "HiReCS" System Flow Cytometry/Cell Sorting at the University of Texas: A Non-Commercial Flow Cytometer

The University of Texas has developed diverse technologies in the field of cytometry allowing them to construct their own flow cytometers. The University of Texas facility consists of a home-built, High-Resolution Cell Sorter (the "HiReCS" system) equipped with two lasers. A dual argon-ion laser/dye laser system (a 5 watt UV-enhanced, water-cooled Spectra/Physics Model 2025 argon-ion laser and a 300mW air-cooled argon-ion laser (Ion Laser Technology Model 5500AWC) and a Spectra/Physics Model 350 dye head unit which is optically split to provide up to three excitation laser beams at two spatial locations. This configuration provides up to 6-color fluorescence analysis capabilities. Complex 3- and 4-color fluorescence experiments are routinely performed.

Unlike most commercial systems, a home-built computer-controlled 3-color compensation system allows for simultaneous correction of three overlapping fluorescent signals. A single-cell sorting module allows for sterile cloning of individual cells and allows sorting of single cells directly into tubes for subsequent PCR enzymatic amplification of DNA or RNA. A number of other home-built, patented devices are incorporated into the system. These include a high-resolution time-of-flight measurement system capable of sizing cells and subcellular organelles as small as 0.3 microns in length or diameter; a multiparameter high-speed "rare-event" analysis system ("HISPEED") capable of analyzing cells at rates in excess of 100,000 cells per second. With this device, analysis of rare cell subpopulations as small as 0.0001 percent with sample rates of more than 10⁹ cells per hour can be achieved.

A home-built sophisticated data acquisition system with Windows™ graphical user interface runs on a 66 MHz 80486 computer and allows for acquisition of 8-parameter, 12-bit listmode data first into RAM, then onto a 200 megabyte disk at rates in excess of 10,000 cells per second. The advanced and unique HiReCS also permits visualization, analysis, and sorting of cells from very complex multi-color fluorescence experiments. Correlation between multiple flow cytometric parameters can be explored through the use of multidimensional principal component/biplot analyses and associated real-time sorting. These analysis and sorting methods can also be used to reduce false-positives. A unique expert system has been developed providing "flexible sorting" strategies to

optimize sample yield or purity. Instrumentation is continually designed or updated in response to particular biological applications.

Partec™ Flow Cytometers (CyFlow® and PAS®)

Partec (<http://www.partec.de>) has developed the CyFlow®, a compact flow cytometer for diverse flow cytometry analysis including the capability to determine absolute volumetric counting. The CyFlow allows to analyze FSC, SSC and up to three fluorescent signals. Since the FloMax® software is PC based, the data acquisition, analysis, and real-time display is performed with any type of PC or laptop within a windows platform. The software allows postacquisition compensation, a very important feature during the analysis.

PAS® is a more sophisticated Compact Flow Cytometer from Partec and it is designed for Multiparameter Flow Cytometry analysis for applications in Immunology,, Microbiology, Biology, Pathology, Hematology including Lymphocyte Subset and Cell biology Analysis

Equipped with an Argon ion laser (488 nm) and an HBO lamp for UV excitation it can determine up to 8 parameters: FSC, SSC, FL1 to FL6 (FITC, PI, DAPI, SR, TR, PE-CY5, AMCA, APC, and many more colors and combinations)

It can detect less than 100 FITC molecules, and determine the true volumetric absolute counting that is very useful for assess CD 34+, CD4/CD8 and other immunolabelled cells without need for reference particles.

The Partec flow cytometers are very robust and compact as well as easy to handle and less demanding in location space, maintenance and service. In addition, the fact that some Partec flow cytometers like the CyFlow® runs on 12 V DC power (e.g. on a car battery) can be important for users in remote places or in field research and for developing countries.

Non-coding RNAs and Epigenetic Regulation

Drivers of Natural Selection

Edited by: KV Morris

c. 220 pp, February 2012

ISBN: 978-1-904455-94-3, \$310/£159

An important and up-to-date overview of the modulation of gene transcription by non-coding RNAs. An essential reference book and a major information resource for those working in the area.

Brucella

Molecular Microbiology and Genomics

Edited by: I López-Goñi, DO' Callaghan

c. 288 pp, February 2012

ISBN: 978-1-904455-93-6, \$310/£159

Highly acclaimed *Brucella* scientists comprehensively review the most important advances in the field. Topics include: genetic diversity, proteomic analysis, transcriptomic analysis, and much more.

Flaviviruses

Edited by: P-Y Shi

c. 410 pp, January 2012

ISBN: 978-1-904455-92-9, \$350/£180

An up-to-date and cutting-edge anthology from the leading experts in the flavivirus field. Essential reading for flavivirus researchers at the graduate level and beyond.

Bacterial Pathogenesis

Molecular and Cellular Mechanisms

Edited by: C Locht, M Simonet

c. 380 pp, January 2012

ISBN: 978-1-904455-91-2, \$350/£180

Distinguished scientists comprehensively describe the most relevant and up-to-date information on pathogenic features across the bacterial world.

Bunyaviridae

Molecular and Cellular Biology

Edited by: A Plyusnin, RM Elliott

c. 250 pp, September 2011

ISBN: 978-1-904455-90-5, \$310/£159

Up-to-date reviews of current research in the field of Bunyaviridae molecular and cellular biology.

Antibacterial Discovery

Answering the Call to Arms

Edited by: AA Miller, PF Miller

c. 480 pp, August 2011

ISBN: 978-1-904455-89-9, \$350/£180

A major reference volume on antibacterial research and how it impacts on public health worldwide, the book is essential reading for everyone working in antibacterial research.

ORDER FROM (UK/Europe): Caister Academic Press, c/o Book Systems Plus, 1st Floor, 8 Hill St., Saffron Walden, Essex, CB10 1JD, UK
Tel: 01799 524458 Fax: 01799 524459 <http://uk.caister.com>

ORDER FROM (USA): Caister Academic Press, c/o ISBS, Inc., 920 NE 58th Avenue, Suite 300, Portland OR 97213-3786, USA
Tel: 503 287-3093 Fax: 503 280-8832 <http://usa.caister.com>

Epigenetics

A Reference Manual

Edited by: JM Craig, NC Wong

c. 450 pp, September 2011

ISBN: 978-1-904455-88-2, \$350/£180

A major new volume of 26 comprehensive and cutting-edge reviews. A key reference manual for everyone involved in epigenetics, DNA methylation, cancer epigenetics and related fields and an essential information resource on the fast-moving field of epigenetics.

Metagenomics

Current Innovations and Future Trends

Edited by: D Marco

c. 310 pp, September 2011

ISBN: 978-1-904455-87-5, \$310/£159

The most innovative and recent advances in theoretical, methodological and applied areas of metagenomics. Aimed at researchers and environmental managers involved in metagenomics, students starting research in this field and teachers interested in the new developments.

Nitrogen Cycling

Molecular Analysis

Edited by: JWB Moir

c. 270 pp, July 2011

ISBN: 978-1-904455-86-8, \$310/£159

With contributions from expert authors from around the world, this excellent book provides comprehensive reviews of current nitrogen cycle research.

Paramyxoviruses

Edited by: SK Samal

c. 550 pp, July 2011

ISBN: 978-1-904455-85-1, \$350/£180

A timely and comprehensive review of current paramyxovirus research. Essential reading.

Helicobacter pylori

Edited by: L Boyanova

c. 316 pp, July 2011

ISBN: 978-1-904455-84-4, \$310/£159

Current knowledge and recent research for microbiologists, clinicians and advanced students working with *Helicobacter* and for those wishing to enter the field.

Bioremediation of Nonmetals

Current Research

Edited by: A-I Koukkou

x + 266 pp, July 2011

ISBN: 978-1-904455-83-7, \$310/£159

An essential reference resource for everyone interested in the bioremediation of organic pollutants.

Lactic Acid Bacteria

Current Progress in Advanced Research

Edited by: K Sonomoto, A Yokota

x + 286 pp, July 2011

ISBN: 978-1-904455-82-0, \$310/£159

Essential reading for every researcher working with LAB, bifidobacteria and probiotics, from the PhD student to the experienced scientist.

Viruses and Interferon

Current Research

Edited by: K Mossman

x + 266 pp, May 2011

ISBN: 978-1-904455-81-3, \$310/£159

A panel of international experts reviews current research on the interaction between viruses and interferon.

Veterinary Parasitology

Edited by: HM Elsheikha, NA Khan

x + 222 pp, May 2011

ISBN: 978-1-904455-79-0 (pbk), \$99/£59

ISBN: 978-1-904455-80-6 (hbk), \$169/£99

"an original review" (*Parasites and Vectors*)

Hepatitis C

Antiviral Drug Discovery and Development

Edited by: S-L Tan, Y He

x + 390 pp, April 2011

ISBN: 978-1-904455-78-3, \$350/£180

A timely and well-compiled review of current HCV research with an emphasis on antiviral drug development.

Streptomyces

Molecular Biology and Biotechnology

Edited by: P Dyson

xii + 258 pp, March 2011

ISBN: 978-1-904455-77-6, \$310/£159

"a must-read" (*David A. Hopwood, John Innes Centre, UK*); "a major reference" (*Doodys*)

Alphaherpesviruses

Molecular Virology

Edited by: SK Weller

x + 448 pp, March 2011

ISBN: 978-1-904455-76-9, \$350/£180

This up-to-date and comprehensive volume distills the most important research in this area providing a timely overview of the field.

FULL DETAILS OF ALL OUR BOOKS AT
WWW.CAISTER.COM

Flow Cytometry Resource Center Flow Cytometry Principle Flow Cytometry Sample Preparation Flow Cytometry Protocol Flow Cytometry Optimization Tips Flow Cytometry troubleshooting Tips. Molecular Biology Resources. Molecular Biology Resource Center Molecular Biology Principle Molecular Biology Protocol Molecular Biology Troubleshooting Tips Molecular Biology FAQ. Boster Troubleshooting eBook Series. Want to have all technical references by your finger tips? Download our FREE eBooks for WB, IHC, ELISA, Flow cytometry and Molecular biology here. These ebooks contain detailed information regarding In Flow Cytometry and Cell Sorting. A. Radbruch, ed. Heidelberg: Springer Berlin, 1992Google Scholar.Â FernÃndez A.B., SÃnchez-Moreiras A.M., de la PeÃ±a T.C. (2001) Flow Cytometry: Principles and Instrumentation. In: Reigosa Roger M.J. (eds) Handbook of Plant Ecophysiology Techniques. Springer, Dordrecht. General Principles. Flow cytometry measures optical and fluorescence characteristics of single cells (or any other particle, including nuclei, microorganisms, chromosome preparations, and latex beads). Physical properties, such as size (represented by forward angle light scatter) and internal complexity (represented by right-angle scatter) can resolve certain cell populations.Â Antibodies against intracellular proteins such as myeloperoxidase and terminal deoxynucleotidyl transferase are also commercially available and permit analysis of an increasing number of intracellular markers. erythrocyte analysis. The use of flow cytometry for the detection and quantification of fetal red cells in maternal blood has increased in recent years.